

# **Earth and Environmental Sciences**

## **Dynamics**

### **New IR microscope and bench installed at BL 1.4**

Martin, M.C., H.-Y. Holman, W.R. McKinney

### **SR-FTIR study of bacteria-water interactions: Acid-base titration and silification experiments**

Yee, N., L.G. Benning, M.J. Tobin, K.O. Konhauser

## **New IR microscope and bench installed at BL1.4**

Michael C. Martin<sup>a</sup>, Hoi-Ying N. Holman<sup>b</sup>, and Wayne R. McKinney<sup>a</sup>

<sup>a</sup> Advanced Light Source Division, Lawrence Berkeley National Laboratory, Berkeley, CA 94720

<sup>b</sup> Center for Environmental Biotechnology, Lawrence Berkeley National Laboratory, Berkeley, CA 94720

### **1. INTRODUCTION**

New infrared spectromicroscopy equipment was purchased for and installed on the ALS infrared beamlines on beam port 1.4. It includes the latest step-scan capable FTIR bench and an infinity corrected infrared microscope which will allow for a number of new sample visualization methods. This equipment was purchased with funding from the DOE Office of Biological and Environmental Research (OBER) with the express purpose to develop biomedical and biological applications of synchrotron-based infrared spectromicroscopy.

### **2. EQUIPMENT**

The new spectromicroscopy equipment includes a Thermo Nicolet Nexus 870 step- and rapid-scan FTIR bench, and a Thermo Spectra-Tech Continuum IR microscope, photographed below. The IR microscope includes two IR detectors, a wide-band MCT and a fast (20 ns) TRS MCT for time-resolved experiments. A fast digitizer (up to 100MHz) compliments the TRS MCT detector. The synchrotron beam coupled into the IR microscope continues to have a diffraction-limited spot size, thereby attaining a 200-fold increase in signal from small (3 – 10 micron) sample spot compared to a conventional thermal IR source. The infinity-corrected microscope optics allow for a number of additional sample visualization accessories which can help the user identify the important location within their sample for micro-IR analysis:

- Visual and IR polarizers
- Dark-field illumination
- DIC (Differential Interference Contrast) optics
- UV Fluorescence



An example of DIC optics enhancing a micrograph of human cheek cells is shown in the photograph to the right. The DIC technique provides a psuedo-3D effect, enhancing the contrast between different thicknesses of an otherwise clear sample. In the image to the right, one can make out the nuclei of the cells (thicker bump near the middle of each cell), whereas this would be difficult using conventional illumination.



This new instrument will aide in user scientific research across many fields. For example, the study of individual living cells, toxic contaminants, bioremediation, protein microcrystals, rhizoids, and forensic evidence will all be enhanced by the additional capabilities of this new SR-FTIR spectromicroscopy system.

### **ACKNOWLEDGEMENTS**

This research was supported by the Office of Science, Office of Biological and Environmental Research, Medical Science Division and the Office of Science, Office of Basic Energy Sciences, Materials Sciences Division, of the U.S. Department of Energy under Contract No. DE-AC03-76SF00098 at Lawrence Berkeley National Laboratory.

Principal investigator: Michael C. Martin, Advanced Light Source Division, LBNL, 510-495-2231, [MCMartin@lbl.gov](mailto:MCMartin@lbl.gov)

# SR-FTIR Study of Bacteria-Water Interactions: Acid-base Titration and Silification Experiments

Nathan Yee<sup>1</sup>, Liane G. Benning<sup>1</sup>, Mark J. Tobin<sup>2</sup>, and Kurt O. Konhauser<sup>1</sup>

<sup>1</sup>School of Earth Sciences, University of Leeds, United Kingdom

<sup>2</sup>Synchrotron Radiation Source, Daresbury Laboratory, United Kingdom

Bacterial surfaces are highly reactive and can strongly affect mass transport in a wide range of geological environments. Bacterial cell walls can adsorb aqueous metal cations, and can act as nucleation surfaces for heterogeneous mineral precipitation. However, the reactions at the bacteria-water interface are poorly understood, primarily due to the difficulty in monitoring such processes *in situ* and *in vivo*. In this study, we use synchrotron radiation-based FTIR to investigate the chemistry of bacterial surfaces with acid/base titration and Si precipitation experiments. The objectives of this research are to identify the reactive surface functional groups and to determine how metal adsorption/precipitation affects the protein and lipid structures of individual living bacterial cells. *In-situ* FTIR experiments were performed on the Infrared beamline 1.4.3 at the Advance Light Source (Lawrence Berkeley National Laboratory), using a Nicolet 760 FTIR bench and a Spectra-Tech Nic-Plan IR microscope. All experiments were performed with flow through fluid cell with BaF<sub>2</sub> and ZnSe windows separated by a 6  $\mu$ m mylar spacer. Acid-base titration and Si precipitation experiments were conducted with both intact cells and isolated bacterial sheaths of *Calothrix* (strain KC97) a filamentous cyanobacteria. Titration experiments with intact bacterial cells show a change in peak position of the carboxylic functional group at  $\sim 1400\text{ cm}^{-1}$  (symmetric vibrational stretching of deprotonated carboxylate groups) from acidic to near-neutral pH (Fig 1). The bacterial silicification experiments indicate a change in peak position at  $\sim 1700 - 1740\text{ cm}^{-1}$ , corresponding to the vibrational C=O stretching of esters groups in the lipid structures of the cell (Fig 2). Previous studies have demonstrated that hydrogen bonding onto carbonyl functional

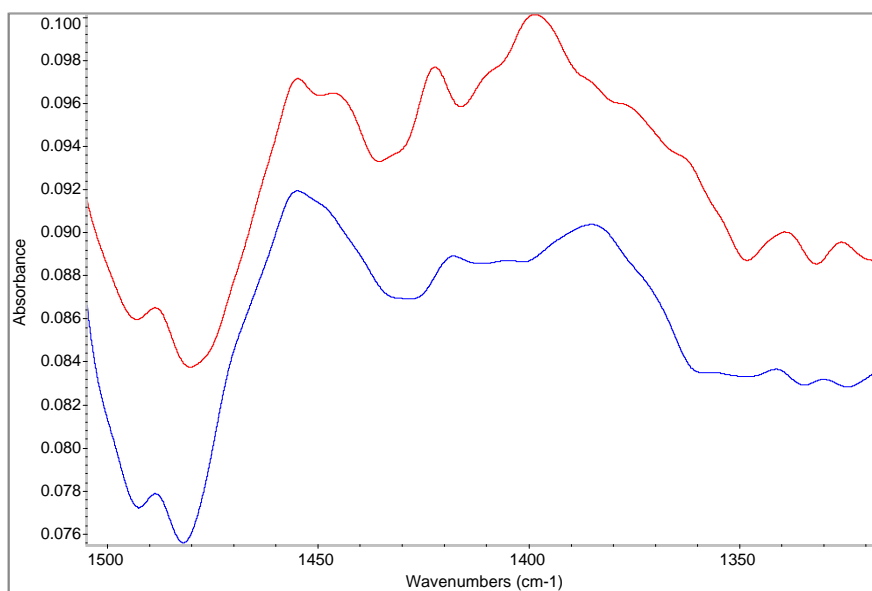


Figure 1. Infrared spectra of an intact *Calothrix* cell in aqueous solution at pH 2.9 and 6.3. A shift in peak position is observed at  $1400\text{ cm}^{-1}$  corresponding to  $\nu_s(\text{COO}^-)$  stretching of deprotonated carboxylate functional groups.

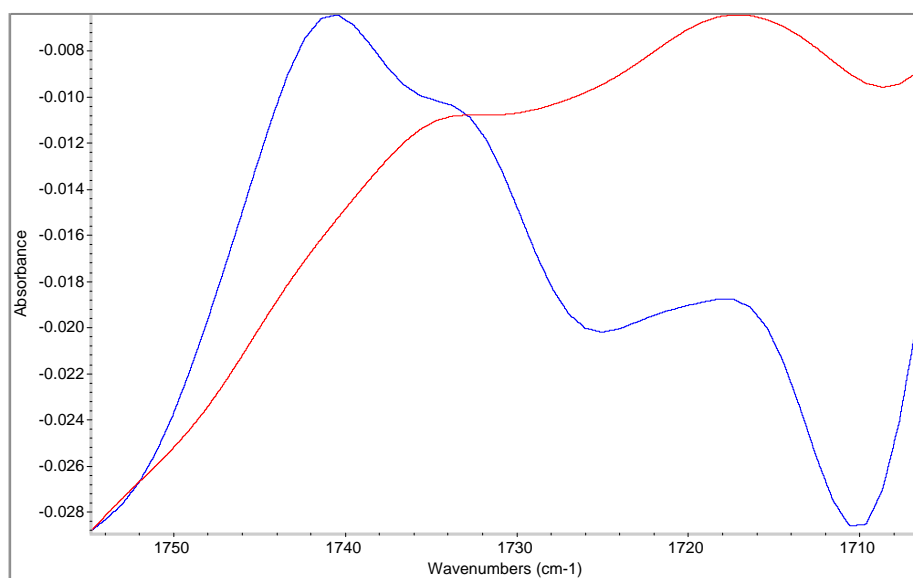


Figure 2. Infrared spectra of an intact *Calothrix* cell before and after silicification. A change in peak position of vibrational C=O stretching is observed.

groups can shift the peak positions in this wavenumber range. Finally, infrared microspectrometry experiments were performed to image the proteins, lipids, and nucleic acids inside intact living cells. Spatial resolution of a few microns was achieved and the chemical distribution of proteins was mapped throughout a *Calothrix* filament (Fig 3). The data indicate that protein molecules have a high concentration within the cell, but a very low concentration on the bacterial surface. These results demonstrate that SR-FTIR can be applied to investigate the functional group chemistry of bacteria in a range of different bacteria-water systems.

This work is an ENVIROSYNC project funded by the Natural Environmental Research Council, UK.

Principal investigator: Liane Benning, University of Leeds, +44-113-233-5220, liane@earth.leeds.ac.uk.

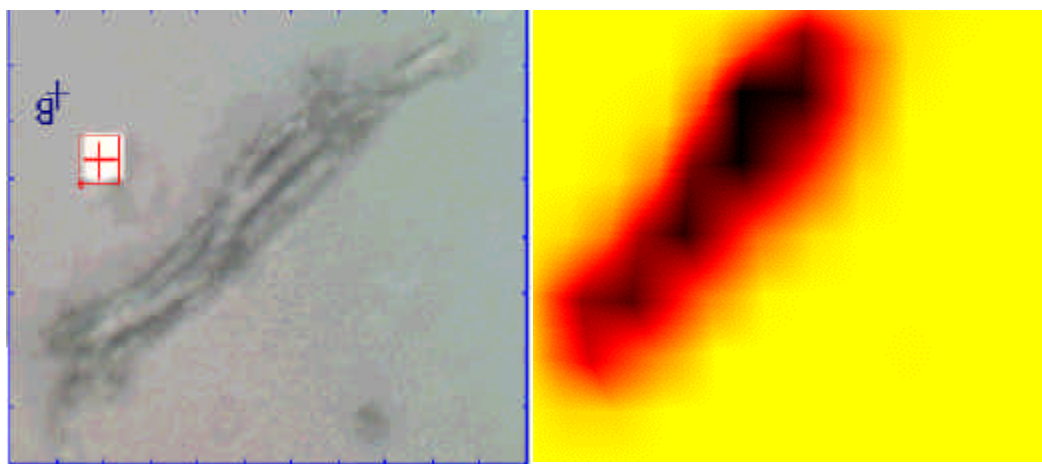


Figure 3. Chemical distribution of the protein characteristic bands amide I and amide II ( $1495\text{--}1800\text{cm}^{-1}$ ) of a *Calothrix* filament. a) Optical image, b) 2-D map of the protein distribution.